Date: 12/21/14 — Page 1/3

Fluorescent Tools for BioMedical Applications

http://www.setabiomedicals.com
e-mail: info@setabiomedicals.com

Product number: K8-1642
Product name: Seta-632-NHS

#### **General Data**

Molecular Mass: 918.02 (protonated form)
Solubility: water, DMF, DMSO

Insoluble: acetone, chloroform, toluene

Storage: Store in absence of light, desiccate and refrigerate

## **Description**

• Highly hydrophilic, amine-reactive label containing one NHS-ester group.

## **Applications**

- · Covalent labeling of proteins, amino-modified DNA and amino-modified oligonucleotides
- Fluorescence intensity and fluorescence polarization-based applications
- Resonance Energy Transfer (RET)
- Flow Cytometry
- Immunofluorescence
- Gene Expression
- · Homogeneous Assays
- · Assessment of protein structure

#### **Advantages**

- Perfectly suited for excitation with the 635 nm diode lasers
- · Sensitive; high extinction coefficients and high quantum yields after covalent attachment to biomolecules
- Quantum yield is highly increased after covalent and non-covalent association with proteins
- pH-insensitive between pH 3 and pH 10
- Good aqueous solubility; this label does not alter the solubility of the bioconjugate
- High photostability; e.g. compared to fluorescein or Cy5<sup>™</sup>
- Low molecular weight Seta dyes do not add substantial mass to the conjugates
- Ideal for non-radioactive labeling of proteins, amino-modified oligonucleotides and amino-modified lipids


http://www.setabiomedicals.com e-mail: info@setabiomedicals.com Product number: K8-1642
Product name: Seta-632-NHS

# **Spectral Data**


Solvent System: phosphate buffer pH 7.4

Sample	Dye-to-protein Ratio	Absorption max. [nm]	Extinction Coefficient [M <sup>-1</sup> cm <sup>-1</sup> ]	Fluorescence max. [nm]	Quantum Yield <sup>1</sup> [%]
Free dye	_	632	280,000	641	6
BSA conjugate 1	1.0	644		654	43
BSA conjugate 2	2.0	643		654	25
BSA conjugate 3	3.0	643		654	17
BSA conjugate 4	4.0	643		654	15
IgG conjugate 1	1.0	637		647	18
IgG conjugate 2	2.0	637		647	13
IgG conjugate 3	3.0	637		647	10
IgG conjugate 4	4.0	637		647	8
IgG conjugate 5	6.0	637		647	7
IgG conjugate 6	8.0	637		647	7


<sup>&</sup>lt;sup>1</sup> Excitation at 600 nm

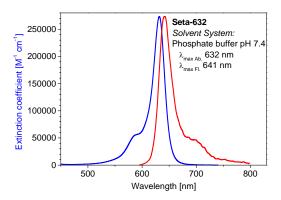

Absorption and emission spectrum of **Seta-632** in phosphate buffer (pH 7.4)

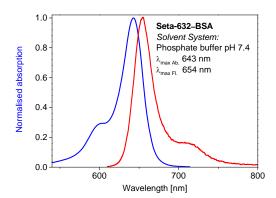

Absorption and emission spectrum of **Seta-632** — **BSA conjugate** in phosphate buffer (pH 7.4, Dye-to-protein ratio 2.0)

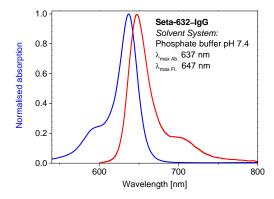


Absorption and emission spectrum of **Seta-632** — **IgG conjugate** in phosphate buffer (pH 7.4, Dye-to-protein ratio 1.0)


Date: 12/21/14 — Page 3/3

Product number: K8-1642
Product name: Seta-632-NHS

http://www.setabiomedicals.com e-mail: info@setabiomedicals.com

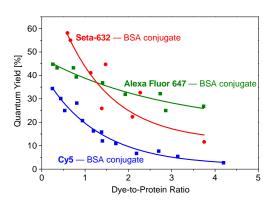

Quantum Yield vs Dye-to-protein Ratio of **Seta-632** — **BSA conjugates** as compared to **Cy5** and **Alexa Fluor 647** in phosphate buffer (pH 7.4)

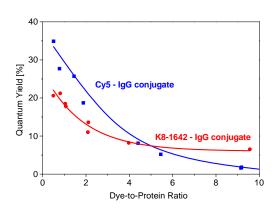

Quantum Yield vs Dye-to-protein Ratio of **Seta-632** — **IgG conjugates** as compared to **Cy5** in phosphate buffer (pH 7.4)